

Stellar M22

Shine on



The future
of skin
treatments
starts now

Shine brightly

There are as many people as there are stars in the sky.
Each one different.
Each one special.

For anyone on the planet, the Stellar M22™
can change their world.

With their skin rejuvenated to the best it can be,
they will feel good, they will feel confident.
They will feel like a star.

Shine on.



With the Stellar M22, the sky is the limit

The Stellar M22 is a powerful modular multi-application platform, inspired by the brightest constellation in the sky.

It enables you to safely and effectively treat different indications in different skin types, age, and gender, without the need for disposables.

The Stellar M22 will empower you to treat more patients and expand your business, creating a distinct competitive advantage for your practice.

The system is fully upgradeable, so as your business grows, the Stellar M22 can adapt and grow with it.



Offer a truly holistic treatment approach



Attract more patients



Benefit from a rapid return on investments

Enable everyone to shine

The Stellar M22 enabling you to treat more patients and expand your business. Discover the infinite possibilities available with the Stellar M22 for these key treatments and many more:

Photofractional™

Sequential treatment of IPL and ResurFX™ for both skin texture and tone, for overall improved skin appearance.

Pigmentation & Tattoos

Non-invasive treatment to reduce unwanted pigmentations and tattoos.

Textural Treatments

Fractional skin rejuvenation solution to improve skin texture.

Vascular Lesions

Non-invasive treatment for Leg Veins, Facial Veins, Broken Capillaries, Rosacea and redness.

Acne

A gentle IPL treatment that addresses both the p.acnes bacteria, as well as the pigmentation and redness of the skin.

One modular system. A galaxy of treatments



One versatile applicator to treat over 20 skin conditions and hair removal

- Treats Acne, Rosacea, Pigmented and Vascular Lesions. Also reduce age spots through photorejuvenation.
- Optimal Pulse Technology (OPT®) controls pulse shape and length.
- Multiple Sequential Pulsing (MSP™) enables a safer procedure while cooling between pulses.
- 9 ExpertFilters™, changeable in seconds.
- 3 SapphireCool™ lightguides for patient comfort through continuous contact cooling.




Needs just one pass for effective skin resurfacing, saving time and protecting patient's skin

- The only true fractional non-ablative technology.
- State-of-the-art CoolScan™ scanner with continuous contact cooling for increased patient comfort.
- Choose from over 600 combinations of area shape, size and density for optimal treatment, whether a single lesion or the whole face.
- Easy visibility of the treatment area through the SapphireCool™ Tip window.



For precise treatment of vascular lesions, and leg veins on all skin types

- Multiple Sequential Pulsing (MSP™) allows cooling between pulses. Reducing risk of skin damage when treating Vascular Lesions: the series of sub-pulses enables deeper penetration and an overall longer pulse duration, which is crucial during treatment.
- 3 SapphireCool™ lightguides, changeable in seconds, for patient comfort through continuous contact cooling.

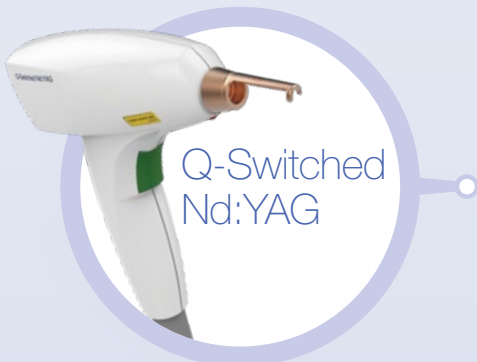


Skin toning by treating pigmented lesions and removing dark tattoos

- Homogenous 'Top Hat' beam profile for enhanced treatment safety and efficacy.
- Single head for all Q-Switched applications.
- 7 spot sizes for customized treatments.

Out of this world technology

The Stellar M22 treats all skin types, with over 30 indications for use. Ability to use multiple cutting-edge technologies in the same treatment session for outstanding clinical results.


OPT® - Optimal Pulse Technology

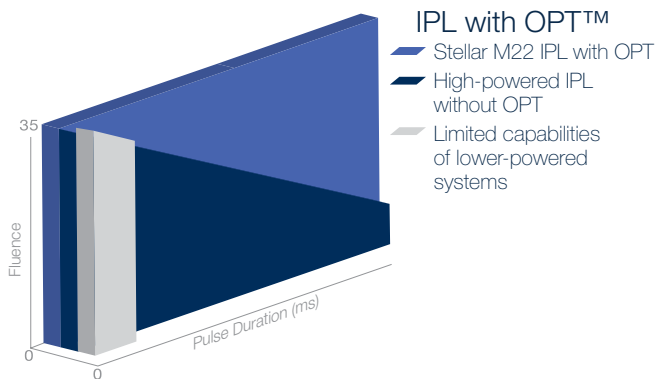

The OPT® advantages:

- High peak power, shorter pulses - optimal for IPL skin treatments using photorejuvenation and the treatment of benign pigmented lesions.
- Homogeneous power over long pulse durations.

MSP™ - Multiple-Sequential Pulsing

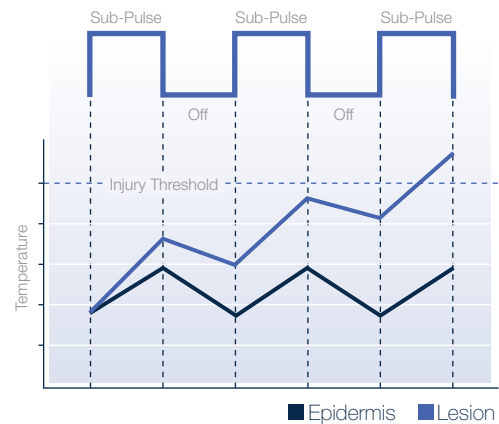

MSP™ Technology:

Multiple-Sequential Pulsing allows for cooling the skin in between pulses, reducing the chance of skin damage during treatments.


ResurFX™ CoolScan™ Shapes

6 different shapes, with sizes and densities changeable from the user interface, allow treatment precision and control



Easy to use and learn

With intuitive user interface, including treatment presets.

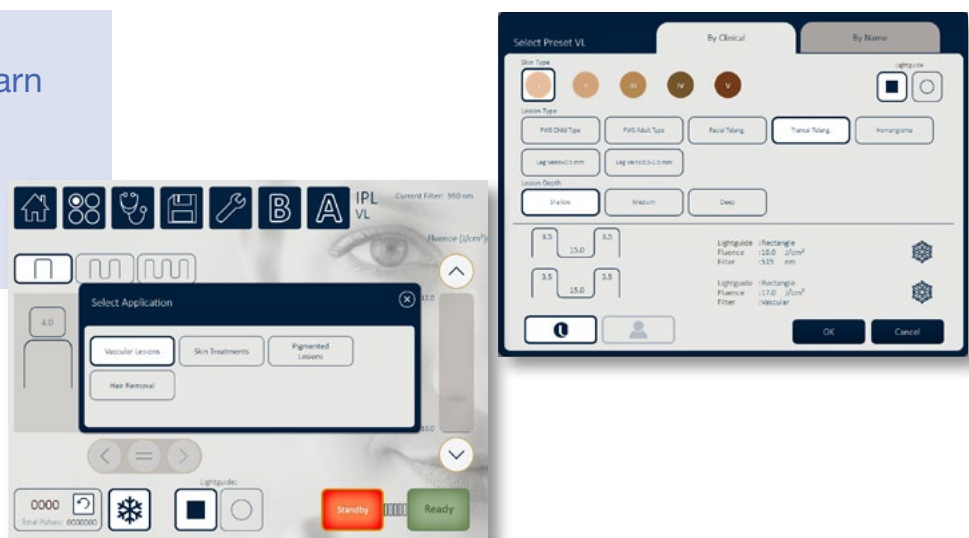

Enable everyone to shine

We're confident in our results, so you can be too

The multi-platform capabilities of the Stellar M22 have been validated in numerous clinical studies and peer-reviewed articles. Practitioners all over the world have continuously reported high levels of patient satisfaction and repeatable results.

Before



After 4 Treatments




Acne treatment with IPL
Courtesy of Li Yuanhong, MD

Before



After

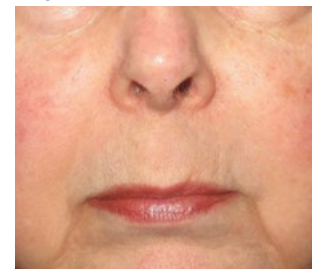

Rosacea treatment with IPL
Courtesy of Gilly Munavalli, MD

Before




After




Skin resurfacing with ResurFX
Courtesy of Matteo Tretti Clementoni, MD

Before

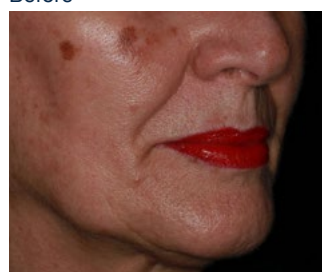

After

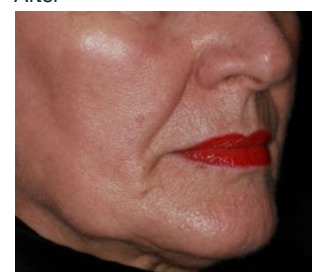

Photofractional™ treatment with IPL & ResurFX
Courtesy of Matteo Tretti Clementoni, MD

Before




After



Leg veins treatment with Nd:YAG
Courtesy of Jose Miguel Montero Garcia, MD

Before

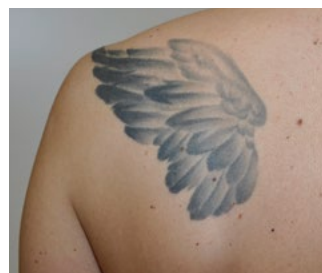

After




Tattoo removal with Nd:YAG Q-Switched
Courtesy of Michael Oeser, MD

Before



After

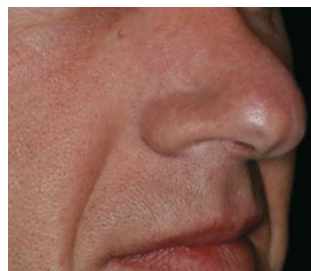

Vascular treatment with IPL
Courtesy of Matteo Tretti Clementoni, MD

*When I think about purchasing a new device, I think about the ROI, and how quickly I'm going to be able to pay for the device... And I was happy to see that **we paid for the device in just few months, which I think is unbelievable!***

Melanie D. Palm, MD, MBA

Your journey with Lumenis has just begun

We want your business to grow with Lumenis. That's why we have put together a support package that begins with the first conversation and will stay with you as your business begins to see the benefits. We provide a support package which ensures that you have all the training you need and receive the most up to date clinical materials. Our direct to consumer marketing activities support the brand, keeping Stellar M22 the treatment of choice when it comes to laser and light skin treatment.



Lumenis Practice Locator

Thousands of patients search for treatment providers every month. Help them find your clinic by registering your practice on the Lumenis Practice Locator. We run campaigns through influencers, social media and viral videos to promote this directory.



Lumenis online

Consumers have many questions about aesthetic. Our patient-facing websites will give them the re-assurance and confidence they need to move forward with the treatments.



Team Training

We install the system and train your staff, so that you are able to assess skin types and conditions, creating the perfect treatment plan for each individual patient. Our expert clinical training means they will handle consultations and treatments with confidence.



Your Most Valuable Resource

Your dedicated representative is an email, text, or phone call away. We want you to feel part of the Lumenis family, so please feel free to contact us any time. We'll always do our best to get back to you within 24 hours.



Brand Your Clinic

As a Lumenis client you will have access to the Lumenis Physician Portal. Here, you can find materials such as brochures and posters. You will also have access to a wealth of clinical and product information.



Extensive Service Network

We are committed to providing high quality, comprehensive service for medical equipment repair and maintenance. Lumenis service provides you with peace of mind, so you can focus on patient care, rather than on the equipment.

Specifications

Stellar IPL™

Spectrum	400-1200 nm
ExpertFilters	Acne (400-600 & 800-1200); Vascular (530-650 & 900-1200); 515 nm, 560 nm; 590 nm; 615 nm; 640 nm; 695 nm; 755 nm;
Fluence	up to 56 J/cm2
Pulse Duration	4-20 ms
Pulse Delay	5-150 ms
Pulse Characteristic	Multiple Sequential Pulsing
Repetition Rate	Up to 1 Hz
Spot Size	35x15 mm ² ; 15x8 mm ² , 6mm
Cooling	Continuous contact cooling

MultiSpot™ Nd:YAG

Wavelength	1064 nm
Fluence	10-225 J/cm ²
Pulse Duration	2-20 ms
Pulse Delay	5-100 ms
Pulse Characteristic	Multiple Sequential Pulsing
Repetition Rate	Up to 1 Hz
Spot Size	6 mm; 2x4 mm ² ; 9 mm
Cooling	Continuous contact cooling

ResurFX™

Wavelength	1565 nm
Pulse Energy	10-70 mJ per micro-beam
Pulse Characteristic	CoolScan™ Scanner
Pattern Shapes	Horizontal line, vertical line, square, horizontal rectangle, vertical rectangle, circle, donut, hexagon
Beam Density	50-500 micro-beams/cm ²
Repetition Rate	0.5-2 Hz
Tip	SapphireCool; Precision
Spot Size	Up to 18 mm
Cooling	Continuous contact cooling

Q-Switched Nd:YAG

Wavelength	1064 nm
Fluence	0.9-14 J/cm ²
Pulse Duration	6-8 ns
Pulse Characteristic	Top-Hat
Repetition Rate	0.5-5.0 Hz
Spot Size	2; 2.5; 3.5; 4; 5; 6 and 8 mm

System specifications

Dimensions [without ResurFX]	51X56X52 cm (21X22X20.5 inches)
Dimensions [with ResurFX]	51X56X65 cm(21X22X25.6 inches)
Weight [without ResurFX]	45 kg / 99.2 lbs
Weight [with ResurFX]	60 kg / 132.3 lbs
Electrical requirements	100 - 240 VAC ±10%, 12A max, 50/60 Hz, single phase, dedicated line

Warnings and risks

As with any laser or IPL, appropriate care must be taken to ensure safe and proper use.

The Stellar M22 operator's manual should be thoroughly reviewed and understood before operating this instrument.

Stellar M22, including the Stellar IPL, Multi-Spot Nd:YAG, Q-Switched Nd:YAG and ResurFX, is contraindicated for patients with active infection, viral, fungal or bacterial diseases. It is also contraindicated for patients with skin disorders or inflammatory skin conditions.

Be advised that IPL/Laser can cause epidermal injury. The risk increases with greater fluence /energy and skin pigmentation.

The use of Stellar M22 could cause redness, swelling, change of pigmentation and scarring. See the system user manual for a complete list of contraindications and risk.

© 2019 Lumenis Ltd., All Rights Reserved. The Lumenis Group of Companies. Lumenis, its logo, Stellar M22, Stellar IPL, OPT, Multiple-Sequential Pulsing, ExpertFilters, SapphireCool lightguides, Photofractional, ResurFX and CoolScan are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice Lumenis is ISO 13485:2012 certified.



Rx ONLY



Lumenis Ltd.

Yokneam Industrial Park
Hakidma Street 6, P.O.B. # 240
Yokneam 2069204, ISRAEL
Tel: +972 4 959 9000

AMERICAS

San Jose, CA, USA
T +1 408 764 3000
+1 877 586 3647
F +1 408 764 3999

JAPAN

Tokyo, Japan
T +81 3 4431 8300
F +81 3 4431 8301

EMEA

Dreieich Dreieichenhain,
Germany
T +49 6103 8335 0
F +49 6103 8335 300

Roma (RM), Italy

T +39 06 90 75 230
F +39 06 90 75 269
Hertfordshire, UK
T +44 20 8736 4110
F +44 20 8736 4119

ASIA / PACIFIC

Beijing, China
T +86 1057 3766 77
F +86 1057 3767 67
Kowloon, Hong Kong
T +852 2174 2800
F +852 2722 5151

Taipei, Taiwan
T +886 2 7719 9098
F: +886 2 7719 9005

Seoul, Korea
T +82 2 6137 9225
F +82 2 6137 9231

Gurgaon, India
T +91 124 4220 795
F +91 124 4854 932

Sydney, Australia
T +61 1800 5863 64

© Lumenis Ltd 2019